

SimpliShade™ BULK FILL



Simplify with One Shade and One Step

Beautiful restorations just got simpler, with SimpliShade™ Bulk Fill. You can skip the capping layer with a one-step placement and reduce inventory requirements with one shade for all 16 VITA® classical shades. Powered by Ternary Adaptive Response Technology (ART³), SimpliShade Bulk Fill is easy to handle with low shrinkage stress.

Kerr™

ONE SHADE IS ALL IT TAKES

Get the adaptability, durability, depth of cure, and simplicity you need with SimpliShade Bulk Fill.



ART³: A NEW WAY TO BULK FILL

Achieve lifelike restorations with the next generation of clinically proven Adaptive Response Technology (ART) found in Harmonize™ and SimpliShade Universal Composite. Powered by Enhanced Adaptive Response Technology, SimpliShade Bulk Fill offers easy handling and low shrinkage.



One-Shade System

Using one shade that blends seamlessly with all 16 VITA® classical shades, SimpliShade Bulk Fill saves you time selecting shades and managing inventory.



One Step

Eliminate the need for a capping layer, complete restorations with a SingleFill™ composite. With its enhanced curing mechanism, SimpliShade Bulk Fill makes it easy to achieve polymerization depths up to 5mm in a single step.




Simplicity with Durability

SimpliShade Bulk Fill offers lower wear rate, higher flexural strength, and comparable shrinkage stress relative to other bulk fill brands on the market, keeping your restorations stronger longer.



Chameleon Effect

Special nanofiller architecture mimics the natural enamel structure and interacts with light for virtually indistinguishable blending.

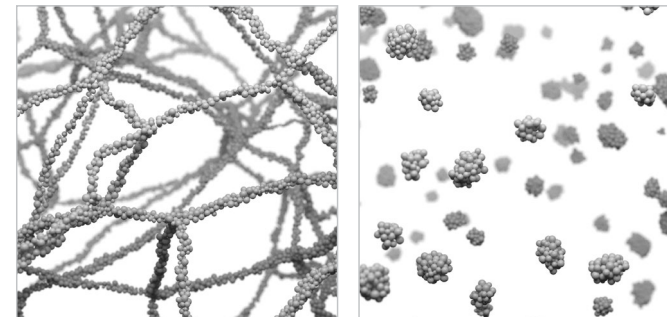


Adaptive Rheological Response Technology for Ease of Sculpting

ART³ uses rheological modifiers that reduce the composite's viscosity as you manipulate it. When you're done shaping, viscosity recovers and anatomy stabilizes before light-curing — reducing the need for resculpting.


Adaptive Color-Blending Response Technology for Beautiful Blending

Achieve exceptional blending with special nanofillers that mimic enamel structures and interact with light in a natural way. The restoration and tooth become virtually indistinguishable.

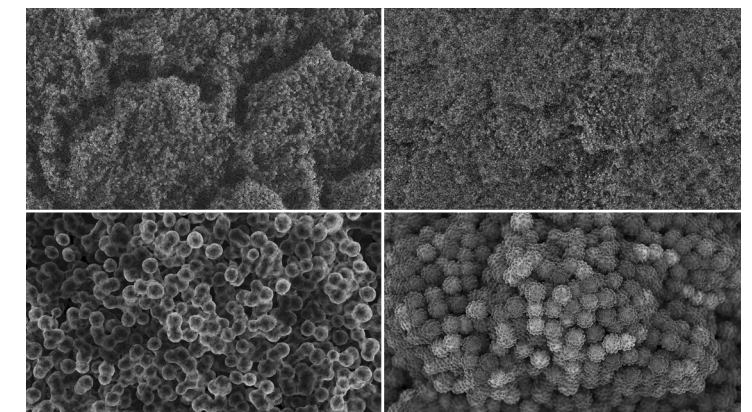

Human Enamel

SimpliShade Bulk Fill

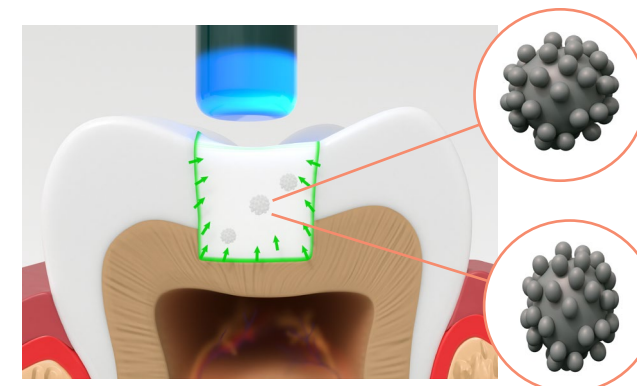


Adaptive Shrinkage Stress Response Technology for Durable Results

Increase confidence in long-lasting restorations. ART³ relieves interfacial stress between the filler and resin matrix during curing. Less stress on the adhesion layer allows for enhanced marginal adaptation.

Flexural Strength

SimpliShade Bulk Fill demonstrates comparable or higher flexural strength compared to Tetric® PowerFill and Admira® Fusion x-tra resulting in better resistance to chipping and fracture over time.

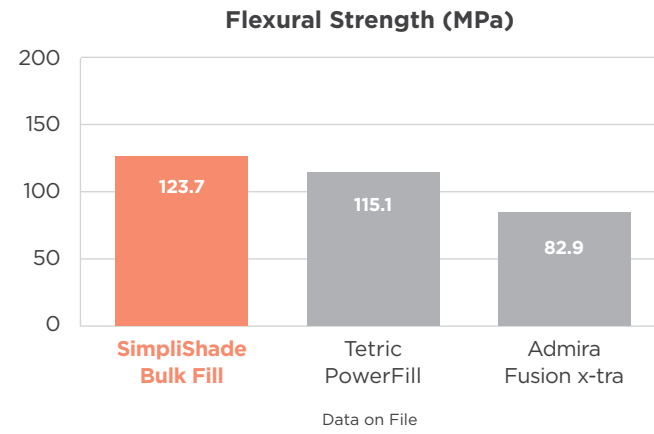


Handling Properties

Adaptive rheological response technology reduces viscosity during manipulation for easier sculpting and anatomy creation.



Shrinkage Stress (at 2 min.)

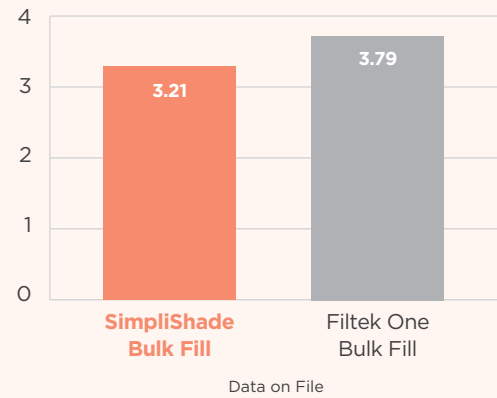


Shrinkage Stress


SimpliShade Bulk Fill has comparable shrinkage stress to Filtek™ One Bulk Fill, resulting in lower marginal degradation and postoperative sensitivity for a long-lasting restoration.

Wear Rate

SimpliShade Bulk Fill has the lowest wear rate relative to Filtek One Bulk Fill, Tetric PowerFill, and Admira Fusion x-tra. Low wear rate results in a restoration that will withstand the daily forces of chewing, brushing, and other occurrences.

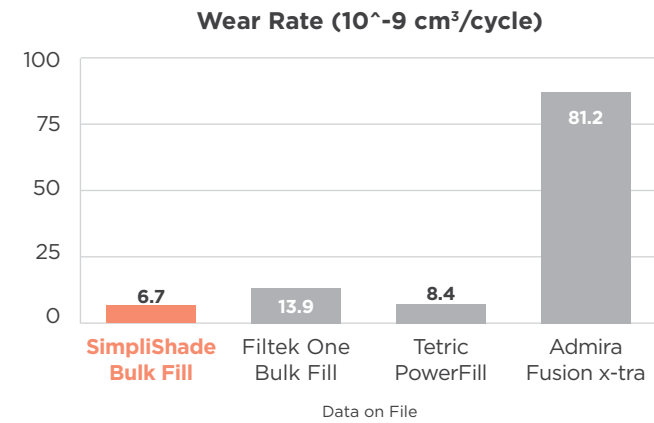

Initial Gloss (GU)

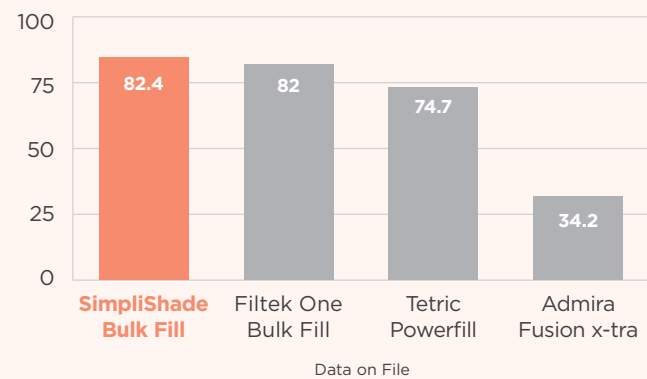

Initial Gloss

SimpliShade Bulk Fill has comparable initial gloss to Filtek One Bulk Fill, and higher initial gloss than Tetric PowerFill and Admira Fusion x-tra, resulting in a more esthetic restoration.

SimpliShade Bulk Fill used here as a capping layer for SimpliShade Bulk Fill Flow.



Ordering Information



SimpliShade Bulk Fill

| Part # | Item Description |
|--------|------------------|
|--------|------------------|

| | |
|-------|-------------------------------------------------|
| 37201 | SimpliShade Bulk Fill Composite Syringe |
| 37202 | SimpliShade Bulk Fill Composite Unidose 20 Pack |

Keep it in the family with our full line of SimpliShade products!



SimpliShade Bulk Fill Flow

| Part # | Item Description |
|--------|------------------|
|--------|------------------|

| | |
|-------|--------------------------------------------|
| 37207 | SimpliShade Bulk Fill Flow Syringe |
| 37208 | SimpliShade Bulk Fill Flow Unidose 20 Pack |



SimpliShade Universal Composite

| Part # | Item Description |
|--------|------------------|
|--------|------------------|

| | |
|-------|-------------------------------------------------------------------|
| 37001 | SimpliShade Universal Composite Syringe, Light |
| 37002 | SimpliShade Universal Composite Syringe, Medium |
| 37003 | SimpliShade Universal Composite Syringe, Dark |
| 37004 | SimpliShade Universal Composite Syringe, Bleach White |
| 37005 | SimpliShade Universal Composite Syringe, Universal Opaque |
| 37006 | SimpliShade Universal Composite Unidose 20 Pack, Light |
| 37007 | SimpliShade Universal Composite Unidose 20 Pack, Medium |
| 37008 | SimpliShade Universal Composite Unidose 10 Pack, Dark |
| 37009 | SimpliShade Universal Composite Unidose 10 Pack, Bleach White |
| 37010 | SimpliShade Universal Composite Unidose 10 Pack, Universal Opaque |



For an ultra high shine finish, use Kerr's line of Spiral Polishers!



One Shade. One Step.
A Better Way to Bulk Fill.

For more information scan the QR Code or visit go.kerrdental.com/simplishadebulkfill